

wiser

LASEREVOLUTION





IT'S TIME TO BELIEVE IN A NEW **LASER.**

Faster treatments, faster healing time, bloodless operations, less anesthesia, less pain: welcome to the perfect dental practice. The support of Wiser diode laser revolutionizes the quality of care for the wellbeing of the patient.

NEW HARMONY BOTH FOR DENTISTS AND PATIENTS.


Fear of dentist is a thing of the past with a Wiser by your side.
Its greatest value is safety.

The patient feels comfortable with the laser and the treatment goes smoothly, without complications. **Wiser is fast, accurate and tireless:** just like the best assistant.

A patient on the way? Keep calm and Wiser will do it all!
The secret is in the technology: Wiser is a surgical laser designed to manage all treatments in **endodontics, periodontics, soft tissue surgery, therapy, biostimulation and whitening.**



Periodontics

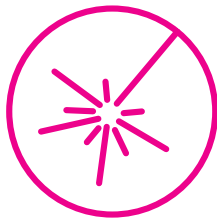

Biostimulation



Surgery



Whitening



wiser



Therapy



Endodontics

THE GREAT REVOLUTION, NOW.

The benefits of Wiser can't be ignored compared to conventional techniques.

Wiser is so gentle that the use of anesthetics is not always necessary.

Think of children and the most difficult patients: this is the real revolution!

Better results, less risk, faster healing, without complications and sutures: this is Tissue Optimized Pulse surgery by Wiser. Greatest comfort with the highest outcomes for the patients.

The revolution is here, now!



MORE VALUE TO YOUR PROFESSIONAL SKILLS.

- Powerful: up to 16W peak of power
- Tissue Optimized Pulsing: specific settings for all types of soft tissue
- 3 selection working modalities: assisted, quick and advanced
- Smart: the Laser Assisted Protocol offers great safety in any procedure
- Intuitive: rapid selection of treatments through simple icons
- Ergonomic: wide touch display
- Portable: it works without cables, with a wireless pedal
- Chargeable: in just 1 minute thanks to the patented SuperCap system
- Autoclavable tips: maximum hygiene and less operating costs
- Stylish: Italian design, innovative and compact

AN EXPERT FOR ANY SPECIALIST.



Performance at its best. The success of the treatment is guaranteed by three different user modes:

- **Laser Assisted Protocol:** a tutorial will guide you step by step with all necessary actions to guarantee the success of any operation. The interface suggests power, time and specific tip to use.
- **Quick mode:** pre-set programs allow fast use with the best golden standard protocols. This mode is designed for faster and simpler use of the laser in most common daily treatments.
- **Advanced mode:** includes a broader range of specific applications in each treatment category. This mode allows the customization of any parameter of the laser emission.

EVOLUTION MEANS KNOWLEDGE.



Do you want to start using laser technology but lack the experience? Would you like to learn more about a specific laser application? Doctor Smile, in collaboration with Master Class Academy, offers courses at various levels held by our skilled laser tutors that include the analysis of different clinical cases and hands on practical sessions. More information and the course agenda are available at masterclassacademy.it



Made in Italy

Wiser – Laserevolution is a CE medical device entirely made in Italy.

Model	Doctor Smile Wiser
Code	LA8D0 001.3
Wavelength	980 nm
Max power	16 W
Pulse mode	Continuous, pulsed, super-pulsed
Frequency	Up to 25 kHz
Size (WxHxD)	161x208x176mm
Weight	1,85 kg
Medical /Laser class	II B / 4

Attestato con Sistema di Certificazione per la Qualità secondo standard



All LAMBDA S.p.A. medical laser are certified with CE trademark in accordance with EEC Directives for medical device: EEC Directives 93/42

A M DP I74.2

Dealer



MORE INFO
wiserlaser.com • doctor-smile.com • info@doctor-smile.com

**DOCTOR SMILE IS
A TRADEMARK OF
LAMBDA SPA**
Via dell'Impresa, 1
36040 BRENDOLA (VI) Italy
T +39 0444.349165
F +39 0444.349954
info@lambdaspa.com
lambdaspa.com